

The DAFO 5 Softy uses a soft foam liner inside a custom-fabricated, UCBL-style foot orthosis to provide comfort and flexible containment of the forefoot and heel. Two toe shelf trimlines are available: the standard option is more flexible, with the plantar surface plastic shell trimmed proximal to the 1st met head, and distal to the 5th, with the Softy liner extending full length. The outer frame can also be extended to full length with the liner, if preferred. Submalleolar trimlines offer a solution for patients who require more control than a shoe insert but don't need the support of an above-ankle orthosis.



## Key Features

**Custom-fabricated from a cast**

**Low profile dorsal trimlines for ease of shoe fit**

**Partial heel and midfoot containment**

**Arch fill and medial heel stabilization**


**Defined plantar surface contours**

**Softy Liner**

provides more comfortable control for difficult-to-correct feet with sensitive bony anatomy


**Easy toe-off distal trimlines**

Plastic outer frame trimmed proximal to 1st met. head, distal to 5th

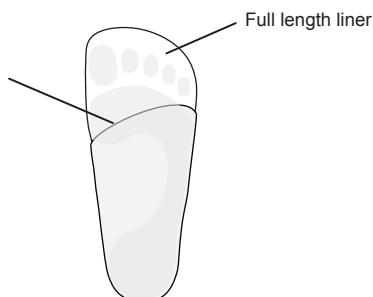


## Designed for patients who:

- Present with moderate pronation and associated gait instability.
- Need more containment than is possible with a lower trimline shoe insert.
- Require extra padding for feet with sensitive bony prominences.

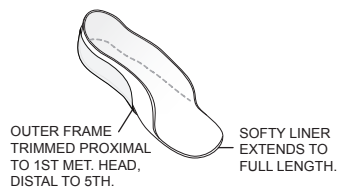

Standard option: More flexible, only liner extends full length

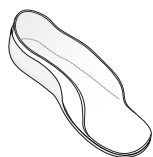

Plastic outer frame and liner extends full length